




- mittelfristig resorbierbares Multifilament
  - hohe Reißkraft
  - minimales Gewebetrauma
  - glatter Gewebedurchzug
  - guter Knotensitz
- *midterm absorbable multifilament*
  - *high tensile strength*
  - *minimized tissue trauma*
  - *smooth tissue passage*
  - *perfect knot security*

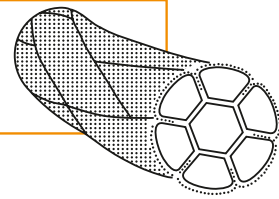

## Eigenschaften:

- hohe Reißkraft
- gutes Handling
- guter Knotensitz
- atraumatischer Gewebedurchzug
- minimale Gewebereaktion

PGA RESORBA® ist ein Polymer der Glykolsäure. Die linear hochmolekulare Polyglykolsäure wird unter Einwirkung eines Katalysators über das Zwischenprodukt Glykolid, einem zyklischen Ester, synthetisiert.

Der Abbau des PGA-Fadenmaterials im Gewebe erfolgt durch Wasseraufnahme in Umkehrung der Synthese.

Die monomere Glykolsäure wird durch den normalen Stoffwechsel enzymatisch in CO<sub>2</sub> und H<sub>2</sub>O gespalten. Fadenmaterial mit einem Anteil von 10% Lactid als Copolymer unterscheidet sich in seinem physikalischen und physiologischen Verhalten nur geringfügig von reinen PGA-Fäden.

Die feinen, präzisionsgeflochtenen Filamente gewähren eine sehr hohe Reißfestigkeit sowie eine ausgezeichnete Geschmeidigkeit.

Die spezielle Beschichtung legt sich fein um die Faserbündel und führt zu einer gezielten Herabsetzung der Oberflächenreibung. Resorbierbare Nahtmaterialien adaptieren das Gewebe während der Heilungsphase und verlieren dabei zunehmend ihre Zugfestigkeit und Reißkraft.

Bei PGA RESORBA® führen die präzisionsgeflochtenen Filamente aus Polyglykolsäure zu einem standardisierten, mittelfristigen Resorptionsverhalten im Gewebe. 18 Tage nach Implantation liegt die Reißkraft von PGA RESORBA® bei 50% der ursprünglichen Reißkraft.\*

Violettes PGA RESORBA® ist mit einem physiologisch unbedenklichen Farbstoff gefärbt.

- ➔ Farbe: violett oder ungefärbt
- ➔ Chemische Bezeichnung: Polyglykolsäure
- ➔ Fadenstärke: USP 8-0 -5 (0.4-7 metric)
- ➔ Konfektionsangebote:
  - Nadel-Faden-Kombinationen
  - in vorgeschrittenen Längen
- ➔ Sterilisation: Ethylenoxid

\* Bei Artikeln, die von Healthium hergestellt werden, liegt die ursprüngliche Reißkraft 21 Tage nach Implantation bei 55%

## Properties:

- high tear resistance
- good handling
- excellent knot security
- atraumatic tissue passage
- minimal tissue reaction

PGA RESORBA® is a polymer of glycolic acid. The linear, high molecular weight polyglycolic acid is synthesised in the presence of a catalyst via the intermediary product glycolide, a cyclic ester.

The breakdown of the PGA suture in the tissue occurs through water

uptake, thereby reversing the synthesis. Regular metabolic processes break down the glycolic acid monomers into CO<sub>2</sub> and H<sub>2</sub>O through enzymatic degradation. The physical and physiological properties of suture material containing 10% lactide as copolymer differ only slightly from pure PGA sutures.

The thin, precision-braided filaments provide very high tear resistance as well as excellent suppleness.

The special coating thinly encapsulates the fibre bundles resulting in a specific reduction of surface friction. Absorbable sutures approximate the tissue during the healing phase; at the same time, they increasingly lose their tensile strength and tear resistance.

PGA RESORBA® contains precision-braided filaments made of polyglycolic acid that result in a standardised, moderate rate of absorption in the tissue. Depending on the thread diameter 50% of the original tensile strength of PGA RESORBA® is lost after 18 days.\*


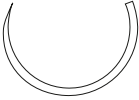




Violet PGA RESORBA® is coloured with a physiologically harmless dye.

- ➔ Colour: violet or undyed
- ➔ Chemical name: polyglycolic acid
- ➔ Thread diameter: USP 8-0 -5 (0.4-7 metric)
- ➔ Types of packaging:
  - needle-thread-combinations
  - precut lengths
- ➔ Sterilization method: ethylene oxide

\* Articles manufactured by Healthium are experience loss of strength down to 55% after 21 days







# PGA RESORBA®

Sterile Nadel-Faden-Kombination / *needle suture combination, sterile*

Nadel / needle	Stärke / gauge		Länge / length (m)	Farbe / colour	Best.-Nr. code no.	Packungsinhalt box content (Stück / pieces)
	USP	metric				
<b>FR = 5/8-Kreis Rundkörper</b>						
<i>FR = 5/8 circle, taper point</i>						
<b>FR 27</b> ● 	2-0	3	0,70	violett / violet	PAN1019 <sup>1)</sup>	24
	0	3.5	0,70	violett / violet	PAN1015 <sup>1)</sup>	24
<b>FR 35</b> ● 	1	4	0,70	violett / violet	PAN1014 <sup>1)</sup>	24
<b>HR, HRS, HRX = 1/2-Kreis Rundkörper</b>						
<i>HR, HRS, HRX = 1/2 circle, taper point</i>						
<b>HR 8</b> ● 	7-0	0.5	0,45	violett / violet	PAN10289 <sup>1)</sup>	12
<b>HR 10</b> ● 	5-0	1	0,70	ungefärbt / undyed	PAN10262 <sup>1)</sup>	24
<b>HR 12</b> ● 	6-0	0.7	0,45	violett / violet	PAN10273 <sup>1)</sup>	24
	5-0	1	0,45	violett / violet	PAN10274 <sup>1)</sup>	24
	6-0	0.7	0,70	violett / violet	PAN1021 <sup>1)</sup>	24
	5-0	1	0,70	violett / violet	PAN1022 <sup>1)</sup>	24
	4-0	1.5	0,70	violett / violet	PAN1023 <sup>1)</sup>	24
	6-0	0.7	0,70	ungefärbt / undyed	PAN10242 <sup>1)</sup>	24
	5-0	1	0,70	ungefärbt / undyed	PAN10265 <sup>1)</sup>	24
	4-0	1.5	0,70	ungefärbt / undyed	PAN10266 <sup>1)</sup>	24
<b>HR 17</b> ● 	5-0	1	0,70	violett / violet	PAN1024 <sup>1)</sup>	24
	4-0	1.5	0,70	violett / violet	PAN1025 <sup>1)</sup>	24
	3-0	2	0,70	violett / violet	PAN1026 <sup>1)</sup>	24
	2-0	3	0,70	violett / violet	PAN1027 <sup>1)</sup>	24
	5-0	1	0,70	ungefärbt / undyed	PAN10248 <sup>1)</sup>	24
	4-0	1.5	0,70	ungefärbt / undyed	PAN10243 <sup>1)</sup>	24
	3-0	2	0,70	ungefärbt / undyed	PAN10244 <sup>1)</sup>	24

# PGA RESORBA®






Sterile Nadel-Faden-Kombination / *needle suture combination, sterile*

Nadel / needle		Stärke / gauge		Länge / length (m)	Farbe / colour	Best.-Nr. code no.	Packungsinhalt box content (Stück / pieces)	
		USP	metric					
<b>HR 17 A</b> ● Abziehnadel / control release needle		5-0	1	4 x 0,45	violett / violet	PA10256 (M)	12	
		4-0	1.5	4 x 0,45	violett / violet	PA10287 (M)	12	
		3-0	2	4 x 0,45	violett / violet	PA10288 (M)	12	
		4-0	1.5	8 x 0,45	violett / violet	PA1028 (M)	12	
		3-0	2	8 x 0,45	violett / violet	PA1029 (M)	12	
		4-0	1.5	4 x 0,70	violett / violet	PA10284 (M)	12	
<b>HR 22</b> ●		4-0	1.5	0,70	violett / violet	PAN10210 <sup>1)</sup>	24	
		3-0	2	0,70	violett / violet	PAN10211 <sup>1)</sup>	24	
		2-0	3	0,70	violett / violet	PAN10212 <sup>1)</sup>	24	
		3-0	2	4 x 0,70	violett / violet	PA10283 (M)	12	
		3-0	2	0,70	ungefärbt / undyed	PAN10254 <sup>1)</sup>	24	
<b>HR 22 A</b> ● Abziehnadel / control release needle		4-0	1.5	4 x 0,45	violett / violet	PA10252 (M)	12	
		3-0	2	4 x 0,45	violett / violet	PA10286 (M)	12	
		2-0	3	4 x 0,45	violett / violet	PA10276 (M)	12	
		4-0	1.5	8 x 0,45	violett / violet	PA10213 (M)	12	
		3-0	2	8 x 0,45	violett / violet	PA10214 (M)	12	
		2-0	3	8 x 0,45	violett / violet	PA10215 (M)	12	
		3-0	2	4 x 0,70	violett / violet	PA10241 (M)	12	
		4-0	1.5	8 x 0,70	violett / violet	PA10285 (M)	12	
<b>HRS 22</b> ●		2-0	3	0,70	violett / violet	PAN10331 <sup>1)</sup>	24	
		0	3.5	0,70	violett / violet	PAN10332 <sup>1)</sup>	24	
		1	4	0,70	violett / violet	PAN10333 <sup>1)</sup>	24	
<b>HRS 22 A</b> ● Abziehnadel / control release needle		0	3.5	4 x 0,70	violett / violet	PA10357 (M)	12	
<b>HR 27</b> ●		4-0	1.5	0,70	violett / violet	PAN10216 <sup>1)</sup>	24	
		3-0	2	0,70	violett / violet	PAN10217 <sup>1)</sup>	24	
		2-0	3	0,70	violett / violet	PAN10218 <sup>1)</sup>	24	
		0	3.5	0,70	violett / violet	PAN10219 <sup>1)</sup>	24	
		4-0	1.5	0,90	violett / violet	PAN10220 <sup>1)</sup>	24	
		3-0	2	0,90	violett / violet	PAN10221 <sup>1)</sup>	24	
		2-0	3	0,90	violett / violet	PAN10222 <sup>1)</sup>	24	
		0	3.5	0,90	violett / violet	PAN10223 <sup>1)</sup>	24	
		3-0	2	0,70	ungefärbt / undyed	PAN10250 <sup>1)</sup>	24	
		2-0	3	0,70	ungefärbt / undyed	PAN10251 <sup>1)</sup>	24	
		0	3.5	0,70	ungefärbt / undyed	PAN10253 <sup>1)</sup>	24	
		2-0	3	6 x 0,70	violett / violet	PA10249 (M)	12	

(M) = Multipack

# PGA RESORBA®








Sterile Nadel-Faden-Kombination / *needle suture combination, sterile*

Nadel / needle	Stärke / gauge USP	metric	Länge / length (m)	Farbe / colour	Best.-Nr. code no.	Packungsinhalt box content (Stück / pieces)	
<b>HR 27 A</b>  Abziehnadel / control release needle	4-0	1.5	4 x 0,45	violett / violet	PA10227 (M)	12	
	3-0	2	4 x 0,45	violett / violet	PA10228 (M)	12	
	2-0	3	4 x 0,45	violett / violet	PA10229 (M)	12	
	4-0	1.5	8 x 0,45	violett / violet	PA10224 (M)	12	
	3-0	2	8 x 0,45	violett / violet	PA10225 (M)	12	
	2-0	3	8 x 0,45	violett / violet	PA10226 (M)	12	
	3-0	2	0,70	violett / violet	PA10280	24	
	4-0	1.5	4 x 0,70	violett / violet	PA10230 (M)	12	
	3-0	2	4 x 0,70	violett / violet	PA10231 (M)	12	
	2-0	3	4 x 0,70	violett / violet	PA10232 (M)	12	
	0	3.5	4 x 0,70	violett / violet	PA10261 (M)	12	
	3-0	2	4 x 0,45	ungefärbt / undyed	PA10263 (M)	12	
	2-0	3	4 x 0,45	ungefärbt / undyed	PA10264 (M)	12	
	<b>HRS 27</b> 	2-0	3	0,70	violett / violet	PAN1030 <sup>1)</sup>	24
0		3.5	0,70	violett / violet	PAN1031 <sup>1)</sup>	24	
1		4	0,70	violett / violet	PAN1032 <sup>1)</sup>	24	
2		5	0,70	violett / violet	PAN10334 <sup>1)</sup>	24	
2-0		3	0,90	violett / violet	PAN10335 <sup>1)</sup>	24	
0		3.5	0,90	violett / violet	PAN1033 <sup>1)</sup>	24	
1		4	0,90	violett / violet	PAN1034 <sup>1)</sup>	24	
2		5	0,90	violett / violet	PAN1035 <sup>1)</sup>	24	
1		4	4 x 0,70	violett / violet	PA10364 (M)	12	
2		5	4 x 0,70	violett / violet	PA10351 (M)	12	
2-0		3	0,70	ungefärbt / undyed	PAN10353 <sup>1)</sup>	24	
<b>HRS 27 A</b>  Abziehnadel / control release needle		2-0	3	4 x 0,45	violett / violet	PA10336 (M)	12
		0	3.5	4 x 0,45	violett / violet	PA10337 (M)	12
	1	4	4 x 0,45	violett / violet	PA10338 (M)	12	
	0	3.5	8 x 0,45	violett / violet	PA10348 (M)	12	
	1	4	8 x 0,45	violett / violet	PA10372 (M)	12	
	0	3.5	4 x 0,70	violett / violet	PA10339 (M)	12	
	1	4	4 x 0,70	violett / violet	PA10340 (M)	12	
	2	5	4 x 0,70	violett / violet	PA10341 (M)	12	
	<b>HRX 27</b> 	1	4	0,90	violett/violet	PAN1041 <sup>1)</sup>	24
<b>HRX 27 A</b>  Abziehnadel / control release needle	1	4	4 x 0,70	violett / violet	PA1045 (M)	12	

(M) = Multipack

# PGA RESORBA®








Sterile Nadel-Faden-Kombination / *needle suture combination, sterile*

Nadel / needle	Stärke / gauge USP	metric	Länge / length (m)	Farbe / colour	Best.-Nr. code no.	Packungsinhalt box content (Stück / pieces)
<b>HR 30</b> ● 	3-0	2	0,70	violett / violet	PAN10245 <sup>1)</sup>	24
	2-0	3	0,70	violett / violet	PAN10233 <sup>1)</sup>	24
	0	3.5	0,70	violett / violet	PAN10234 <sup>1)</sup>	24
<b>HR 30 A</b> ● Abziehnadel / control release needle 	0	3.5	4 x 0,70	violett / violet	PA10270 (M)	12
<b>HRX 30</b> ● 	1	4	0,90	violett / violet	PAN1042 <sup>1)</sup>	24
<b>HR 35</b> ● 	3-0	2	0,70	violett / violet	PAN10246 <sup>1)</sup>	24
	2-0	3	0,70	violett / violet	PAN10235 <sup>1)</sup>	24
	0	3.5	0,70	violett / violet	PAN10236 <sup>1)</sup>	24
	1	4	0,70	violett / violet	PAN10237 <sup>1)</sup>	24
	3-0	2	0,90	violett / violet	PAN10247 <sup>1)</sup>	24
	2-0	3	0,90	violett / violet	PAN10238 <sup>1)</sup>	24
	0	3.5	0,90	violett / violet	PAN10239 <sup>1)</sup>	24
	1	4	0,90	violett / violet	PAN10240 <sup>1)</sup>	24
<b>HR 35 A</b> ● Abziehnadel / control release needle 	3-0	2	4 x 0,45	violett/violet	PA10277 (M)	12
	2-0	3	8 x 0,70	violett/violet	PA10269 (M)	12
	0	3.5	4 x 0,70	violett/violet	PA10267 (M)	12
<b>HRS 35</b> ● 	2-0	3	0,70	violett / violet	PAN1037 <sup>1)</sup>	24
	0	3.5	0,70	violett / violet	PAN1038 <sup>1)</sup>	24
	1	4	0,70	violett / violet	PAN1039 <sup>1)</sup>	24
	3-0	2	0,90	violett / violet	PAN10371 <sup>1)</sup>	24
	2-0	3	0,90	violett / violet	PAN10310 <sup>1)</sup>	24
	0	3.5	0,90	violett / violet	PAN10311 <sup>1)</sup>	24
	1	4	0,90	violett / violet	PAN10312 <sup>1)</sup>	24
	2	5	0,90	violett / violet	PAN10313 <sup>1)</sup>	24
	2-0	3	0,90	ungefärbt / undyed	PAN10349 <sup>1)</sup>	24
<b>HRS 35 A</b> ● Abziehnadel / control release needle 	0	3.5	4 x 0,45	violett / violet	PA10314 (M)	12
	1	4	4 x 0,45	violett / violet	PA10315 (M)	12
	2	5	4 x 0,45	violett / violet	PA10316 (M)	12
	0	3.5	8 x 0,50	violett / violet	PA10324 (M)	12
	1	4	4 x 0,70	violett / violet	PA10343 (M)	12
	2	5	4 x 0,70	violett / violet	PA10344 (M)	12
	0	3.5	4 x 0,45	ungefärbt / undyed	PA10360 (M)	12
	1	4	4 x 0,45	ungefärbt / undyed	PA10361 (M)	12

(M) = Multipack

# PGA RESORBA®

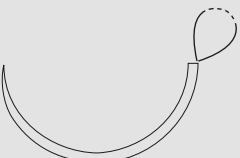

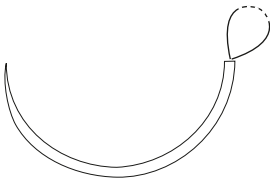


Sterile Nadel-Faden-Kombination / *needle suture combination, sterile*

Nadel / needle	Stärke / gauge USP	Stärke / gauge metric	Länge / length (m)	Farbe / colour	Best.-Nr. code no.	Packungsinhalt box content (Stück / pieces)
<b>HR 40</b> ● 	2-0	3	0,90	violett / violet	PAN10282 <sup>1)</sup>	24
<b>HRS 40</b> ● 	0	3.5	0,70	violett / violet	PAN10317 <sup>1)</sup>	24
	1	4	0,70	violett / violet	PAN10318 <sup>1)</sup>	24
	0	3.5	0,90	violett / violet	PAN10319 <sup>1)</sup>	24
	1	4	0,90	violett / violet	PAN10320 <sup>1)</sup>	24
	2	5	0,90	violett / violet	PAN10321 <sup>1)</sup>	24
<b>HRS 40</b> ● Schlinge / loop 	1	4	1,50	violett / violet	PAN10322 <sup>1)</sup>	24
<b>HRS 40 A</b> ● 	4	1	4 x 0,45	violett / violet	PA10323 (M)	12
	2	5	4 x 0,45	violett / violet	PA10345 (M)	12
	2-0	3	8 x 0,45	violett / violet	PA10365 (M)	12
	1	4	8 x 0,45	violett / violet	PA10373 (M)	12
	2	5	4 x 0,70	violett / violet	PA10375 (M)	12
Abziehnadel / control release needle						
<b>HRS 48</b> ● 	2-0	3	0,70	violett / violet	PAN10368 <sup>1)</sup>	24
	0	3.5	0,70	violett / violet	PAN10377 <sup>1)</sup>	24
	0	3.5	0,90	violett / violet	PAN10325 <sup>1)</sup>	24
	1	4	0,90	violett / violet	PAN10326 <sup>1)</sup>	24
	2	5	0,90	violett / violet	PAN10327 <sup>1)</sup>	24
	2	5	1,20	violett / violet	PAN10354 <sup>1)</sup>	24
	2	5	0,90	ungefärbt / undyed	PAN10374 <sup>1)</sup>	24
<b>HRS 48</b> ● Schlinge / loop 	1	4	1,50	violett / violet	PAN10352 <sup>1)</sup>	24
<b>HRS 48 A</b> ● 	1	4	4 x 0,45	violett / violet	PA10346 (M)	12
	2	5	4 x 0,45	violett / violet	PA10347 (M)	12
	2	5	8 x 0,45	violett / violet	PA10369 (M)	12
Abziehnadel / control release needle						

(M) = Multipack

# PGA RESORBA®

Sterile Nadel-Faden-Kombination / *needle suture combination, sterile*

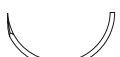
Nadel / <i>needle</i>	Stärke / <i>gauge</i>		Länge / <i>length</i> (m)	Farbe / <i>colour</i>	Best.-Nr. <i>code no.</i>	Packungsinhalt <i>box content</i> (Stück / <i>pieces</i> )
	USP	metric				
<b>HR 50</b> 	1	4	0,70	violett / <i>violet</i>	PAN10278 <sup>1)</sup>	24
	0	3.5	0,90	violett / <i>violet</i>	PAN10279 <sup>1)</sup>	24
<b>HRT, HRTX = 1/2-Kreis schneidende Rundkörper (Diamond-Spitze)</b> <i>HRT, HRTX = 1/2 circle, taper cutting</i>						
<b>HRT 18</b> 	4-0	1.5	0,70	violett / <i>violet</i>	PAN1060 <sup>1)</sup>	24
	3-0	2	0,70	violett / <i>violet</i>	PAN1061 <sup>1)</sup>	24
	5-0	1	0,70	ungefärbt / <i>undyed</i>	PAN10640 <sup>1)</sup>	24
	4-0	1.5	0,70	ungefärbt / <i>undyed</i>	PAN10621 <sup>1)</sup>	24
	3-0	2	0,70	ungefärbt / <i>undyed</i>	PAN10613 <sup>1)</sup>	24
<b>HRT 18 A</b>  <i>Abziehnadel / control release needle</i>	2-0	3	6 x 0,45	violett / <i>violet</i>	PA10620 (M)	12
<b>HRT 22</b> 	3-0	2	0,70	violett / <i>violet</i>	PAN10632 <sup>1)</sup>	24
	3-0	2	0,70	ungefärbt / <i>undyed</i>	PAN10625 <sup>1)</sup>	24
<b>HRT 22 A</b>  <i>Abziehnadel / control release needle</i>	2-0	3	4 x 0,45	violett / <i>violet</i>	PA1068 (M)	12

(M) = Multipack



# PGA RESORBA®



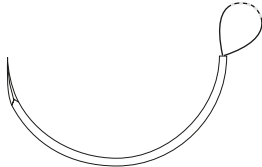


Sterile Nadel-Faden-Kombination / *needle suture combination, sterile*

Nadel / <i>needle</i>	Stärke / <i>gauge</i>		Länge / <i>length</i> (m)	Farbe / <i>colour</i>	Best.-Nr. <i>code no.</i>	Packungsinhalt <i>box content</i> (Stück / <i>pieces</i> )
	USP	metric				
<b>HRT 25</b> 	3-0	2	0,70	violett / <i>violet</i>	PAN1064 <sup>1)</sup>	24
	2-0	3	0,70	violett / <i>violet</i>	PAN1065 <sup>1)</sup>	24
	0	3.5	0,70	violett / <i>violet</i>	PAN1066 <sup>1)</sup>	24
	3-0	2	0,90	violett / <i>violet</i>	PAN10641 <sup>1)</sup>	24
	2-0	3	0,90	violett / <i>violet</i>	PAN10642 <sup>1)</sup>	24
	3-0	2	0,70	ungefärbt / <i>undyed</i>	PAN10622 <sup>1)</sup>	24
	2-0	3	0,70	ungefärbt / <i>undyed</i>	PAN10628 <sup>1)</sup>	24
<b>HRT 25 A</b>  Abziehnadel / <i>control release needle</i>	0	3,5	4 x 0,70	violett / <i>violet</i>	PA10614 (M)	12
	3-0	2	4 x 0,45	ungefärbt / <i>undyed</i>	PA10630 (M)	12
	2-0	3	4 x 0,45	ungefärbt / <i>undyed</i>	PA10631 (M)	12
<b>HRT 35</b> 	1	4	0,70	violett / <i>violet</i>	PAN10634 <sup>1)</sup>	24
	2-0	3	0,90	violett / <i>violet</i>	PAN10610 <sup>1)</sup>	24
	0	3,5	0,90	violett / <i>violet</i>	PAN10611 <sup>1)</sup>	24
	1	4	0,90	violett / <i>violet</i>	PAN10612 <sup>1)</sup>	24
	2	5	0,90	violett / <i>violet</i>	PAN10624 <sup>1)</sup>	24
	2-0	3	0,90	ungefärbt / <i>undyed</i>	PAN10623 <sup>1)</sup>	24
<b>HRT 35 A</b>  Abziehnadel / <i>control release needle</i>	1	4	4 x 0,45	violett / <i>violet</i>	PA10619 (M)	12
<b>HRT 40</b> 	0	3,5	0,70	ungefärbt / <i>undyed</i>	PAN10636 <sup>1)</sup>	24
<b>HRT 45</b>  Schlinge / <i>loop</i>	1	4	1,50	violett / <i>violet</i>	PAN10627 <sup>1)</sup>	24

(M) = Multipack

# PGA RESORBA®








Sterile Nadel-Faden-Kombination / *needle suture combination, sterile*

Nadel / needle	Stärke / gauge		Länge / length (m)	Farbe / colour	Best.-Nr. code no.	Packungsinhalt box content (Stück / pieces)
	USP	metric				
<b>HRTX 50</b> 	1	4	0,90	violett / violet	PAN10617 <sup>1)</sup>	24
<b>ART = Asymptotische schneidende Rundkörper</b> <i>ART = asymptotic, taper cutting</i>						
<b>ART 25</b> 	4-0	1,5	0,45	violett / violet	PAN1281 <sup>1)</sup>	12
<b>DR, DRF = 3/8-Kreis Rundkörper</b> <i>DR, DRF = 3/8 circle, taper point</i>						
<b>DR 12 A</b> 	3-0	2	4 x 0,45	violett / violet	PA1072 (M)	12
<b>DRT = 3/8-Kreis schneidende Rundkörper</b> <i>DRT = 3/8 circle, taper cutting</i>						
<b>DRT 12</b> 	6-0	0,7	0,45	violett / violet	PAN1081 <sup>1)</sup>	24
	5-0	1	0,45	violett / violet	PAN1082 <sup>1)</sup>	24
<b>DRT 18</b> 	4-0	1,5	0,45	ungefärbt / undyed	PAN1083 <sup>1)</sup>	24

(M) = Multipack







# PGA RESORBA®

Sterile Nadel-Faden-Kombination / *needle suture combination, sterile*

Nadel / <i>needle</i>	Stärke / <i>gauge</i>		Länge / <i>length</i> (m)	Farbe / <i>colour</i>	Best.-Nr. <i>code no.</i>	Packungsinhalt <i>box content</i> (Stück / <i>pieces</i> )
	USP	metric				
<b>HS, HSX = 1/2-Kreis außen schneidend</b> <b>HSM = mit PREMIUM-Schliff</b> <i>HS, HSX = 1/2 circle, reverse cutting</i> <i>HSM = PREMIUM cut</i>						
<b>HSM 13</b>  PREMIUM-Schliff / <i>PREMIUM cut</i>	5-0	1	0,45	ungefärbt / <i>undyed</i>	PAN1120 <sup>1)</sup>	24
<b>HS 18</b> 	4-0	1.5	0,70	violett / <i>violet</i>	PAN11112 <sup>1)</sup>	24
	3-0	2	0,70	violett / <i>violet</i>	PAN11113 <sup>1)</sup>	24
	2-0	3	0,70	violett / <i>violet</i>	PAN11114 <sup>1)</sup>	24
<b>HS 22</b> 	4-0	1.5	0,70	violett / <i>violet</i>	PAN11119 <sup>1)</sup>	24
	3-0	2	0,70	violett / <i>violet</i>	PAN11117 <sup>1)</sup>	24
	2-0	3	0,70	violett / <i>violet</i>	PAN11118 <sup>1)</sup>	24
	3-0	2	0,70	ungefärbt / <i>undyed</i>	PAN11119 <sup>1)</sup>	24
<b>HSX 22</b> 	2-0	3	0,70	violett / <i>violet</i>	PAN11117 <sup>1)</sup>	24
<b>HS 35</b> 	2-0	3	0,70	violett / <i>violet</i>	PAN11110 <sup>1)</sup>	24
	0	3.5	0,70	violett / <i>violet</i>	PAN11111 <sup>1)</sup>	24
	1	4	0,70	violett / <i>violet</i>	PAN11112 <sup>1)</sup>	24
<b>HS 40</b> 	1	4	0,90	violett / <i>violet</i>	PAN11113 <sup>1)</sup>	24
	2	5	0,90	violett / <i>violet</i>	PAN11114 <sup>1)</sup>	24
<b>HS 50</b> 	1	4	0,70	violett / <i>violet</i>	PAN11116 <sup>1)</sup>	24





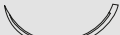

# PGA RESORBA®

Sterile Nadel-Faden-Kombination / *needle suture combination, sterile*

Nadel / needle	Stärke / gauge USP	Stärke / gauge metric	Länge / length (m)	Farbe / colour	Best.-Nr. code no.	Packungsinhalt box content (Stück / pieces)
<b>DS = 3/8-Kreis außen schneidend</b> <b>DSM = mit PREMIUM-Schliff</b> <i>DS = 3/8 circle, reverse cutting</i> <i>DSM = PREMIUM cut</i>						
<b>DSM 11</b>  PREMIUM-Schliff / PREMIUM cut	6-0	0,7	0,45	ungefärbt / undyed	PAN1158 <sup>1)</sup>	24
<b>DS 13</b> 	4-0	1,5	0,70	ungefärbt / undyed	PAN11422 <sup>1)</sup>	24
<b>DSM 13</b>  PREMIUM-Schliff / PREMIUM cut	5-0	1	0,45	violett / violet	PAN1159 <sup>1)</sup>	24
	6-0	0,7	0,70	violett / violet	PAN11512 <sup>1)</sup>	24
	5-0	1	0,70	violett / violet	PAN11511 <sup>1)</sup>	24
	6-0	0,7	0,45	ungefärbt / undyed	PAN1150 <sup>1)</sup>	24
	5-0	1	0,45	ungefärbt / undyed	PAN1151 <sup>1)</sup>	24
	4-0	1,5	0,45	ungefärbt / undyed	PAN1152 <sup>1)</sup>	24
<b>DS 18 F</b> 	4-0	1,5	0,70	ungefärbt / undyed	PAN11414 <sup>1)</sup>	24
<b>DS 18</b> 	5-0	1	0,45	violett / violet	PAN11417 <sup>1)</sup>	24
	4-0	1,5	0,45	violett / violet	PAN11418 <sup>1)</sup>	24
	4-0	1,5	0,70	violett / violet	PAN11423 <sup>1)</sup>	24
	3-0	2	0,70	violett / violet	PAN11416 <sup>1)</sup>	24
	5-0	1	0,45	ungefärbt / undyed	PAN1140 <sup>1)</sup>	24
	4-0	1,5	0,45	ungefärbt / undyed	PAN1141 <sup>1)</sup>	24
	3-0	2	0,45	ungefärbt / undyed	PAN1142 <sup>1)</sup>	24
	4-0	1,5	0,70	ungefärbt / undyed	PAN11424 <sup>1)</sup>	24
	3-0	2	0,70	ungefärbt / undyed	PAN11413 <sup>1)</sup>	24
<b>DSM 18</b>  PREMIUM-Schliff / PREMIUM cut	5-0	1	0,70	violett / violet	PAN11521 <sup>1)</sup>	24
	4-0	1,5	0,70	violett / violet	PAN11510 <sup>1)</sup>	24
	4-0	1,5	0,45	ungefärbt / undyed	PAN1154 <sup>1)</sup>	24
	3-0	2	0,45	ungefärbt / undyed	PAN1155 <sup>1)</sup>	24

# PGA RESORBA®















Sterile Nadel-Faden-Kombination / *needle suture combination, sterile*

Nadel / <i>needle</i>	Stärke / <i>gauge</i> USP	metric	Länge / <i>length</i> (m)	Farbe / <i>colour</i>	Best.-Nr. <i>code no.</i>	Packungsinhalt <i>box content</i> (Stück / <i>pieces</i> )
<b>DS 24</b> ▼ 	2-0	3	0,70	violett / <i>violet</i>	PAN11412 <sup>1)</sup>	24
	4-0	1.5	0,70	ungefärbt / <i>undyed</i>	PAN1143 <sup>1)</sup>	24
	3-0	2	0,70	ungefärbt / <i>undyed</i>	PAN1145 <sup>1)</sup>	24
	2-0	3	0,70	ungefärbt / <i>undyed</i>	PAN1146 <sup>1)</sup>	24
<b>DSM 24</b> ▽ PREMIUM-Schliff / <i>PREMIUM cut</i> 	3-0	2	0,45	ungefärbt / <i>undyed</i>	PAN1156 <sup>1)</sup>	24
	2-0	3	0,45	ungefärbt / <i>undyed</i>	PAN1157 <sup>1)</sup>	24
	3-0	2	0,70	ungefärbt / <i>undyed</i>	PAN11519 <sup>1)</sup>	24
<b>DS 30</b> ▼ 	0	3.5	0,70	violett / <i>violet</i>	PAN11411 <sup>1)</sup>	24
	2-0	3	0,70	ungefärbt / <i>undyed</i>	PAN1149 <sup>1)</sup>	24
	3-0	2	0,70	ungefärbt / <i>undyed</i>	PAN11421 <sup>1)</sup>	24
	0	3.5	0,70	ungefärbt / <i>undyed</i>	PAN1148 <sup>1)</sup>	24
<b>DS 40</b> ▼ 	2-0	3	0,70	ungefärbt / <i>undyed</i>	PAN11420 <sup>1)</sup>	24
<b>GS = Gerade schneidend</b> <i>GS= straight, cutting</i>						
<b>GS 12</b> ▼ 	4-0	1.5	0,70	ungefärbt / <i>undyed</i>	PAN1171 <sup>1)</sup>	24
<b>HSPM, HSPM F = 1/2-Kreis SPATULA, mit PREMIUM-Schliff</b> <i>HSPM, HSPM F = 1/2 circle, spatulated, PREMIUM cut</i>						
<b>2 x HSPM 7</b> ▽ PREMIUM-Schliff / <i>PREMIUM cut</i> 	8-0	0.4	0,30	violett / <i>violet</i>	PAN1195 <sup>1)</sup>	12

\*eventuell mit verlängerter Lieferzeit verbunden / *possibly associated with an extended delivery time*

# PGA RESORBA®

Sterile Nadel-Faden-Kombination / *needle suture combination, sterile*

Nadel / needle	Stärke / gauge		Länge / length (m)	Farbe / colour	Best.-Nr. code no.	Packungsinhalt box content (Stück / pieces)	
	USP	metric					
<b>DSPM = 3/8-Kreis SPATULA, mit PREMIUM-Schliff</b>							
<i>DSPM = 3/8 circle, spatulated, PREMIUM cut</i>							
<b>DSPM 6</b> 		8-0	0.4	0,30	violett / violet	PAN1204 <sup>1)</sup>	12
		7-0	0.5	0,45	violett / violet	PAN1205 <sup>1)</sup>	24
PREMIUM-Schliff / PREMIUM cut							
<b>2 x DSPM 6</b> 		8-0	0.4	0,45	violett / violet	PAN1207 <sup>1)</sup>	12
		7-0	0.5	0,45	violett / violet	PAN1208 <sup>1)</sup>	24
PREMIUM-Schliff / PREMIUM cut							
<b>2 x DSPM 7</b> 		7-0	0.5	0,45	violett / violet	PAN1201 <sup>1)</sup>	24
		PREMIUM-Schliff / PREMIUM cut					
<b>2 x DSPM 8</b> 		6-0	0.7	0,25	ungefärbt / undyed	PAN12014 <sup>1)</sup>	24
		PREMIUM-Schliff / PREMIUM cut					
<b>DSPM 9</b> 		6-0	0.7	0,45	violett / violet	PAN1202 <sup>1)</sup>	24
		PREMIUM-Schliff / PREMIUM cut					
<b>VSPM, VSPM F = 1/4-Kreis SPATULA, mit PREMIUM-Schliff</b>							
<i>VSPM, VSPM F = 1/4 circle, spatulated, PREMIUM cut</i>							
<b>2 x VSPM 8</b> 		5-0	1	0,45	violett / violet	PAN12112 <sup>1)</sup>	24
		6-0	0.7	0,45	violett / violet	PAN12113 <sup>1)</sup>	24
PREMIUM-Schliff / PREMIUM cut							
<b>2 x VSPM 8 F</b> 		6-0	0.7	0,45	violett / violet	PAN1210 <sup>1)</sup>	24
		5-0	1	0,45	violett / violet	PAN1211 <sup>1)</sup>	24
PREMIUM-Schliff / PREMIUM cut							

# PGA RESORBA®

Einzelfadenpackungen, steril in vorgeschrittenen Längen / *Sterile pre-cut suture, non-needed*

Stärke / <i>gauge</i> USP	Länge / <i>length</i> (m)	Farbe / <i>colour</i>	Best.-Nr. <i>code no.</i>	Packungsinhalt <i>box content</i> (Stück / <i>pieces</i> )	
					metric
4-0	1.5	3 x 0,45	violett / <i>violet</i>	PSN333 <sup>1)</sup>	24
3-0	2	3 x 0,45	violett / <i>violet</i>	PSN433 <sup>1)</sup>	24
2-0	3	3 x 0,45	violett / <i>violet</i>	PSN533 <sup>1)</sup>	24
0	3.5	3 x 0,45	violett / <i>violet</i>	PSN633 <sup>1)</sup>	24
1	4	3 x 0,45	violett / <i>violet</i>	PSN733 <sup>1)</sup>	24
2	5	3 x 0,45	violett / <i>violet</i>	PSN833 <sup>1)</sup>	24
4-0	1.5	6 x 0,45	violett / <i>violet</i>	PSN365 <sup>1)</sup>	24
3-0	2	6 x 0,45	violett / <i>violet</i>	PSN465 <sup>1)</sup>	24
2-0	3	6 x 0,45	violett / <i>violet</i>	PSN565 <sup>1)</sup>	24
0	3.5	6 x 0,45	violett / <i>violet</i>	PSN665 <sup>1)</sup>	24
1	4	6 x 0,45	violett / <i>violet</i>	PSN765 <sup>1)</sup>	24
2	5	6 x 0,45	violett / <i>violet</i>	PSN865 <sup>1)</sup>	24
3-0	2	12 x 0,45	violett / <i>violet</i>	PSN4125 <sup>1)</sup>	24
2-0	3	12 x 0,45	violett / <i>violet</i>	PSN5125 <sup>1)</sup>	24
0	3.5	12 x 0,45	violett / <i>violet</i>	PSN6125 <sup>1)</sup>	24
1	4	12 x 0,45	violett / <i>violet</i>	PSN7125 <sup>1)</sup>	24
2	5	12 x 0,45	violett / <i>violet</i>	PSN8125 <sup>1)</sup>	24
3-0	2	2 x 0,70	violett / <i>violet</i>	PSN420 <sup>1)</sup>	24
2-0	3	2 x 0,70	violett / <i>violet</i>	PSN520 <sup>1)</sup>	24
0	3.5	2 x 0,70	violett / <i>violet</i>	PSN620 <sup>1)</sup>	24
1	4	2 x 0,70	violett / <i>violet</i>	PSN720 <sup>1)</sup>	24
2	5	2 x 0,70	violett / <i>violet</i>	PSN820 <sup>1)</sup>	24
4-0	1.5	5 x 0,70	violett / <i>violet</i>	PSN350 <sup>1)</sup>	24
3-0	2	5 x 0,70	violett / <i>violet</i>	PSN450 <sup>1)</sup>	24
2-0	3	5 x 0,70	violett / <i>violet</i>	PSN550 <sup>1)</sup>	24
0	3.5	5 x 0,70	violett / <i>violet</i>	PSN650 <sup>1)</sup>	24
1	4	5 x 0,70	violett / <i>violet</i>	PSN750 <sup>1)</sup>	24
2	5	5 x 0,70	violett / <i>violet</i>	PSN850 <sup>1)</sup>	24

# PGA RESORBA®

Einzelfadenpackungen, steril in vorgeschrittenen Längen / *Sterile pre-cut suture, non-needed*

Stärke / <i>gauge</i> USP	<i>metric</i>	Länge / <i>length</i> (m)	Farbe / <i>colour</i>	Best.-Nr. <i>code no.</i>	Packungsinhalt <i>box content</i> (Stück / <i>pieces</i> )
2-0	3	1 x 1,50	violett / <i>violet</i>	PSN515 <sup>1)</sup>	24
0	3.5	1 x 1,50	violett / <i>violet</i>	PSN615 <sup>1)</sup>	24
1	4	1 x 1,50	violett / <i>violet</i>	PSN715 <sup>1)</sup>	24
2	5	1 x 1,50	violett / <i>violet</i>	PSN815 <sup>1)</sup>	24
4-0	1.5	3 x 0,45	ungefärbt / <i>undyed</i>	PUN333 <sup>1)</sup>	24
3-0	2	3 x 0,45	ungefärbt / <i>undyed</i>	PUN433 <sup>1)</sup>	24
2-0	3	3 x 0,45	ungefärbt / <i>undyed</i>	PUN533 <sup>1)</sup>	24
4-0	1.5	6 x 0,45	ungefärbt / <i>undyed</i>	PUN365 <sup>1)</sup>	24
3-0	2	6 x 0,45	ungefärbt / <i>undyed</i>	PUN465 <sup>1)</sup>	24
2-0	3	6 x 0,45	ungefärbt / <i>undyed</i>	PUN565 <sup>1)</sup>	24
2-0	3	5 x 0,70	ungefärbt / <i>undyed</i>	PUN550 <sup>1)</sup>	24

<sup>1)</sup> Produkt hergestellt durch / *product manufactured by:*

 **Healthium Medtech Limited** · #472-D, 13th Cross, 4th Phase  
Peenya Industrial Area, Bangalore, Karnataka-560058, India  
Email: [care@healthiummedtech.com](mailto:care@healthiummedtech.com) · Mfg. Lic. No.: MFG/MD/2021/000367